SEIDIGYN®

GEL LUBRICANTE ÍNTIMO

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SEIDIGYN®

INTIMATE LUBRICANT GEL

Relieves vaginal dryness favouring sexual intercourse.
External use.

Seidigyn® vaginal lubricant gel is a watery base non-hormonal intimate lubricant gel for external use formulated to help alleviate the symptoms of vaginal dryness, thus favoring sexual intercourse.

Seidigyn® is an odourless, colourless gel with a similar consistency of physiological vaginal secretions, which hydrates and lubricates the genital area reducing friction during sexual intercourse and without compromising condom's integrity.

Vaginal dryness (vulvar and other related disturbances) affects a large number of women. It is usually accompanied by characteristic symptoms such as irritation, vaginal itching and painful intercourse (dyspareunia).

These symptoms may appear:

- In any circumstances causing alterations of estrogenic levels (hormone) in women: pregnancy, post-partum, breastfeeding, contraceptive intake, menopause, chemo or radiotherapy, vaginal surgery and medicines.
- In cases of local hypersensitivity: stress, soap or deodorants allergies or those produced by synthetic or tight underwear.

Seidigyn® intimate lubricant gel is used for the relief of the above mentioned symptoms, which are associated to atrophy and genital dryness. It is respectful with vaginal mucosa and maintains its pH at physiological balance.

It does not contain neither hormonal nor chemical aggressive ingredients.

The soft texture of Seidigyn® forms a lubricant film over the genital surface, providing a comfortable, dryness-relieving feeling when applied over the perineum (vaginal opening, external genital area) and over the vulva.

The lubricant effect of Seidigyn® favours sexual intercourse and it is compatible with the use of condoms.

Seidigyn® vaginal lubricant does not leave stains after vaginal application or over the external genitalia.

Posology

Spread a small amount before each sexual intercourse. Application may be repeated if necessary.

Ingredients

Hydroxyethylcellulose 2%, Glycerol, Propylene glycol, Xanthan gum, Phenoxyethanol, Benzoic acid, Dehydroacetic acid, Inulin, Purified water and Sodium hydroxide.

Contraindications

The product is contraindicated in case of known hypersensitivity to any of its components. Stop using the product in case of hypersensitivity reactions.

Side effects

Not described.

If you suffer any possible side effect while using this product, consult your doctor or pharmacist. Adverse events can be notified through Agencia Española de Medicamentos y Productos Sanitarios portal <https://notificaps.aemps.es/enviotelematico/notificaps/notifica/inicio.do> and sending an email to farmacovigilancia@lab-seid.com. Any serious incident that has occurred in relation to the device should be reported to Seid, S.A. sending an email to farmacovigilancia@lab-seid.com and to the competent authority Agencia Española de Medicamentos y Productos Sanitarios portal <https://notificaps.aemps.es/enviotelematico/notificaps/notifica/inicio.do>

Conservation instructions

- Do not use the product after the expiry date printed on the package.
- Keep the tube closed and away from close contact with heat sources Do not eat.
- Avoid any contact with eyes. In case of accidental contact wash thoroughly with water.
- Keep the product out of the sight and the reach of children.
- Use the product only according to the instructions
- Do not use if the package is damaged.
- Pregnant women and nursing mothers can use the product as recommended by the doctor or pharmacist.
- **Seidigyn® vaginal lubricant** is not a contraceptive and does not contain spermicides. Not prevent pregnancy or sexually transmitted diseases.
- Take care of the environment. Do not throw the packaging in the

Pack size


Aluminium tube containing 30 g of gel.

How to use

Wash your hands before and after using it.

In a Seidigyn® pack you will find an aluminium tube containing 30 g of gel.

1. Unscrew the tube cap turn it down and use the sharp point to perforate the tube's seal.
2. Press the tube so as to release a small amount of gel on your fingertip.
3. Spread the gel over the external intimate area - or over the condom, once it has been put in place, in order to make sexual intercourse more pleasant.
4. The product can be repeatedly applied before each continuous sexual intercourse.

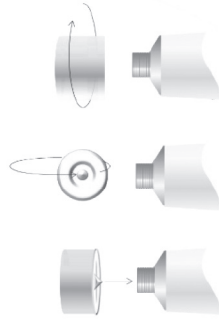Overdosage

There are no risks specifically associated with the overdosage of the product. In case of overdosage remove the excess with a tissue.

Date of last revision: February 2023

Seid, S.A.
Ctra. de Sabadell a Granollers, km 15.
08185 Lliçà de Vall - Barcelona (Spain).

Do not use if package is damaged and consult instructions for use.

Consult instructions for use.

Batch code.

Use-by date.

Unique device identifier

Seidylube 30 g

0051